




## THE CASE FOR AEDS

### Survival increases with early intervention

Research shows that the probability of survival increases dramatically when CPR is performed and when an AED is applied before an ambulance arrives.<sup>1</sup>

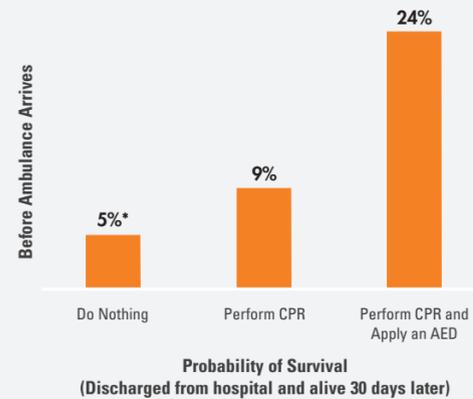

\*US only

### Where is the AED?

Too often, the answer is, "We don't have one." Sadly, the same research that demonstrated a nearly fivefold increase in survivability (from 5% to 24%) when an AED is used, also showed that an AED is available only 2% of the time.

# POWERHEART® G5



#### Intellisense™ Defibrillation Pads

Non-polarized pads can be placed in either location on the chest for simple, fast placement.

#### RescueCoach®

User-paced voice and text prompts guide rescuers through each critical step of the rescue, including CPR instruction.

#### USB Data Download

90 minutes of rescue data available via USB download.

#### Dual Language

Change from primary to secondary language with the push of a button, without delaying or interrupting the rescue.



#### Automatic Self-Testing

Conducts regular, automatic self-tests on pads, battery and internal circuitry to ensure your AED is always Rescue Ready®.

#### Intellisense™ CPR Feedback

Real-time voice and text corrective prompts for CPR quality.

Whether you are an experienced rescuer or a first-time responder, you will appreciate how the Powerheart G5 AED combines ease of use, reliability, and advanced technologies to help victims survive sudden cardiac arrest.

The Powerheart G5 AED with Intellisense CPR provides real-time, Guideline-driven feedback to ensure delivery of high-quality compressions. One-button dual language functionality and comprehensive self-tests give rescuers the confidence to act.

## ZOLL MEDICAL CORPORATION

An Asahi Kasei Group Company | 269 Mill Road | Chelmsford, MA 01824 | 978-421-9655 | 800-804-4356 | www.zoll.com  
 Copyright © 2020 ZOLL Medical Corporation. All rights reserved. Designed for Unexpected Heroes, Rescue Ready, ZOLL, Powerheart, Intellisense and RescueCoach are trademarks or registered trademarks of ZOLL Medical Corporation in the United States and/or other countries. All other trademarks are the property of their respective owners.

<sup>1</sup>Weisfeldt ML, et al. *J Am Coll Cardiol.* 2010;55(16):1713-20.

<sup>2</sup>AHA 2015 Guidelines: pg. S402

<sup>3</sup>AHA CPR Quality Consensus Statement. *Circulation.* 2013; 128; pg. 420

<sup>4</sup>Chung MK, et al. *J Am Coll Cardiol.* 2010 Jul 13;56:194-203. Ong MEH, et al. *JAMA.* 2006;295:22:2629-37

Printed in U.S.A. MCN PP 2010 0516

<sup>5</sup>Kleinman et al. Part 5: adult basic life support and cardiopulmonary resuscitation quality: 2015 American Heart Association Guidelines Update for Cardiopulmonary Resuscitation and Emergency Cardiovascular Care. *Circulation.* 2015;132(suppl 2):S414-S435

<sup>6</sup>*Annals Emergency Medicine; Analysis of Automated External Defibrillator Device Failures, Reported to the Food and Drug Administration, Volume 59, no. 2, February 2012*



## POWERHEART® G5



# DESIGNED FOR UNEXPECTED HEROES™

When someone is the victim of sudden cardiac arrest, their survival depends on high-quality CPR and a rescuer with an AED.

**When there's a call for help, will you be ready?**



# ARE YOU PREPARED TO HELP IN A HEARTBEAT?

It happens in a split second: A person collapses – the victim of sudden cardiac arrest – and the clock starts ticking. You call for help and medical personnel are on their way. But precious seconds are racing by. The victim’s best chance for survival depends on immediate CPR and a rescuer with an automated external defibrillator (AED).



CPR plus early defibrillation can **more than double the rate of survival** from out-of-hospital cardiac arrest.  
-2015 AHA Guidelines<sup>2</sup>

## THE POWER TO HELP

The latest guidelines from resuscitation councils worldwide are clear: Successful defibrillation must be supported by high-quality CPR.

The AED’s first heart analysis is critical. If it calls for no shock, only high-quality CPR can lead to a shockable rhythm on the next heart analysis. High-quality CPR provides a struggling heart with oxygenated blood needed to help restore a normal rhythm.

“...**monitoring of CPR quality is arguably one of the most significant advances in resuscitation practice in the past 20 years and should be incorporated into every resuscitation and every professional rescuer program.**”  
- American Heart Association CPR Quality Consensus Statement<sup>3</sup>

## DESIGNED FOR UNEXPECTED HEROES™

The survival rate from sudden cardiac arrest (SCA) is less than 6% worldwide.<sup>4</sup> According to the American Heart Association (AHA), high-quality CPR, along with early defibrillation, can help improve SCA survival rates.<sup>5</sup>

The Powerheart® G5 AED is the only AED to combine real-time CPR feedback with dualLanguage functionality and user-paced rescue prompting.



“**I was not expecting the AED to walk me through doing compressions. I was sure thankful for that ‘press, press, press’.**”  
Mark Tucker, Rescuer  
CM Tucker Lumber in Pageland, SC



**Real-Time CPR Feedback**  
The Powerheart G5 AED uses Intellisense CPR (ICPR) technology to monitor CPR chest compressions, providing rescuers with distinct corrective voice and text prompts: Press Slower, Press Faster, Press Softer, and Press Harder and Fully Release.



**Rescue Ready®**  
The Powerheart G5 AED performs comprehensive self-testing to ensure it is ready to rescue. Studies have shown that approximately 46% of AED failures are due to pad, pad connector, or battery power problems.<sup>6</sup> The Powerheart AED self-tests detect and help prevent these issues.



**User-Paced Rescue Prompts**  
RescueCoach® prompts guide users step-by-step through the rescue process. The AED recognizes actions taken and ensures that the rescuer completes every critical task before moving on to the next.



**Dual-Language Capabilities**  
With a touch of a button, CPR Instruction, ICPR corrections, and AED prompts can quickly switch from the primary language to a pre-programmed alternate language without delaying the rescue.

