

3870M Large Capacity Manual Sterilizer



3870M

Item No.	Chamber Dimensions Ø x L in(mm)	Volume	Overall Dimensions DxWxH in(mm)	Tray Dimensions DxWxH in(mm)	No. of Trays	Standard Unwrapped Cycle Time	Voltage(V) Freq.(Hz)*	Power(W) Current(A)	Shipping Weight lbs(kg)
3870M with stand	15"x30" (384x760)	22gal/85L	34.5"x26"x20.6" (875x665x525)	26.6"x11.3"x1" (675x286x25) 26.6"x13.8"x1" (675x350x25)	2	Cold: 32 mins Hot: 23 mins	230V 50/60	3000W 13.6A	278 (126) including stand

NOTE: cycle times are approximate and may vary depending on instrument loads. *Model 3870EA requires a voltage between 220V and 235V for proper operation. When applicable a buck/boost transformer (0.5KVA) may be required.

Features

- Automatic shut off at the end of both the sterilization and dry cycles
- A long life electro polished chamber and door
- Double safety locking device prevents door from opening while chamber is pressurized
- Drain valve is located on the front, allowing for quick and easy draining of water reservoir
- Dual safety thermostat to protect against overheating
- Supplied with a 3 shelf stand (30"Wx36"Dx34"H)
- 1 year parts and labor warranty and 10 year chamber warranty

Control Panel



3870EA Large Capacity Automatic Sterilizer



3870EA

Item No.	Chamber Dimensions Ø x L in(mm)	Volume	Overall Dimensions DxWxH in(mm)	Tray Dimensions DxWxH in(mm)	No. of Trays	Standard Unwrapped Cycle Time	Voltage(V) Freq.(Hz)*	Power(W) Current(A)	Shipping Weight lbs(kg)
3870EA with stand	15"x30" (384x760)	22gal/85L	34.5"x26"x20.6" (875x665x525)	26.6"x11.3"x1" (675x286x25) 26.6"x13.8"x1" (675x350x25)	2	Cold: 32 mins Hot: 23 mins	230V 50/60	3000W 13.6A	278 (126) Including Stand

NOTE: cycle times are approximate and may vary depending on instrument loads. *Model 3870EA requires a voltage between 220V and 235V for proper operation. When applicable a buck/boost transformer (0.5KVA) may be required.

Features

- Closed door active drying system to maintain sterility and ensure efficient drying of packs and pouches
- A long life electro polished chamber and door
- New brighter and easier to read panel
- Program parameters can be customized and stored
- Supplied with a 3 shelf stand (30"Wx36"Dx34"H)
- 2 year parts and labor warranty and 10 year chamber warranty

Control Panel



1. Program Selection
2. Detailed Digital Display
3. Programmable Keys
4. System Indicators


Optional Internal Printer



Documents: date, time, temperature, pressure.



Clean & Simple™ Ultrasonic Cleaners and Tablets

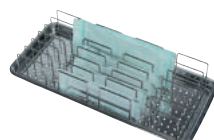

Instrument Pouch Rack



Steam Distiller



Chamber Brite Autoclave Cleaner



Tuttnauer USA Co., Ltd.

25 Power Drive, Hauppauge, NY 11788 Tel: +800 624 5836,
+631 737 4850, Fax: +631 737 0720
E-mail: info@tuttnauerUSA.com, www.tuttnauerUSA.com

